



doctorVet+ therapy laser



EVERYTHING YOU NEED.
EVERYTHING YOU WANT.

3
Keystroke
Protocol
Selection

WiFi
Connectivity

4
Treatment
Head
Options

7" Color
Touch
Screen

Continuous
Wave to
25,000 Hz
635 808 915nm
520nm
(Aiming Beam)

16 W peak power

Battery
Powered

Flexible
2 m Fiber with
Finger Switch



FULLY INTEGRATED TREATMENT MODALITY.

The introduction of **doctorVet Plus** takes Class IV laser therapy in a new direction. Protocols are improved, treatments are enhanced, and User Interface is customizable. The biggest difference is the focus on work flow, connectivity, and DICOM medical records. **doctorVet Plus** has revolutionized the operator, pet owner, case management experience. All the features you need to treat animals effectively, and communicate efficiently with their owners, can be found in the new **doctorVet Plus** Therapy Laser.

MODULAR, EXPANDABLE, SELECTABLE, PLATFORM.

doctorVet Plus is a EWP – Expandable Workflow Platform

Lambda/doctorVet understands practice needs are different around the world. The Expandable Workflow Platform (EWP) allows practices to customize doctorVet Plus to meet their needs. Practices may select:

- Simple Patient Treatments
- Client and Patient Database with treatment history
- Customized Patient Treatment Reports sent to pet owners
- Connectivity to MWL (Modality Worklist) from Practice Management Software
- Connectivity to PACS for organized patient treatment record

Designing doctorVet Plus with an Expandable Workflow Platform shows Lambda understands practices require full customization to meet their clinical, management, and client communication needs.




WORKFLOW & CONNECTIVITY INTEGRATION

doctorVet Plus integrates connectivity providing features and benefits unique to laser therapy. Built-in WiFi and Bluetooth allows for:

Feature	Benefit
Automatic Software Updates	Up To Date Software Features (Just Like Your Smartphone)
Personalized Patient Treatment Reports	Improved Client Communications and Compliance
DICOM connection to PACS	Improved Record Keeping & Legal Medical Record
Modality Work List (MWL) from practice management software	Improved Workflow

Each feature can be customized by individual practices.



DOSIMETRY AND PROTOCOL DEVELOPMENT FOR OPTIMAL PATIENT CARE.

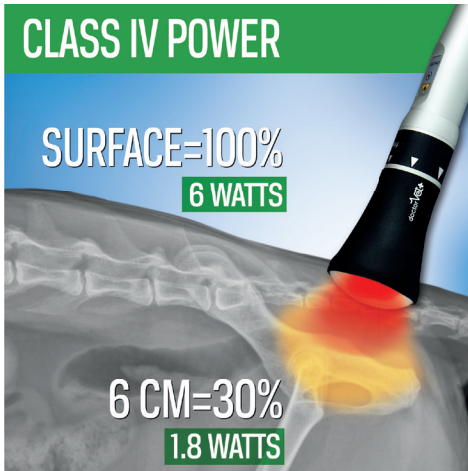


Optimal patient care results from incorporating Power, Pulsing, and the science of Dosimetry. While it is true that Continuous Wave (CW) delivers the greatest dosage in the fastest time, studies prove that variable pulsing also provides enhanced stimulation to specific tissues.

Pulsing above 10,000 Hz is also shown to provide increased anti-bacterial and anti-microbial effects. By combining 16W peak power and 25,000 Hz pulsing capabilities with scientifically supported protocols, **doctorVet Plus** consistently ensures treatment efficacy.


STIMULATING WAVELENGTHS.

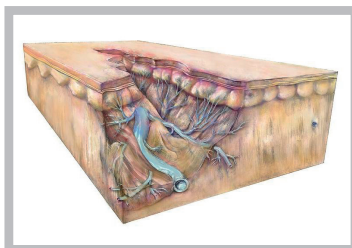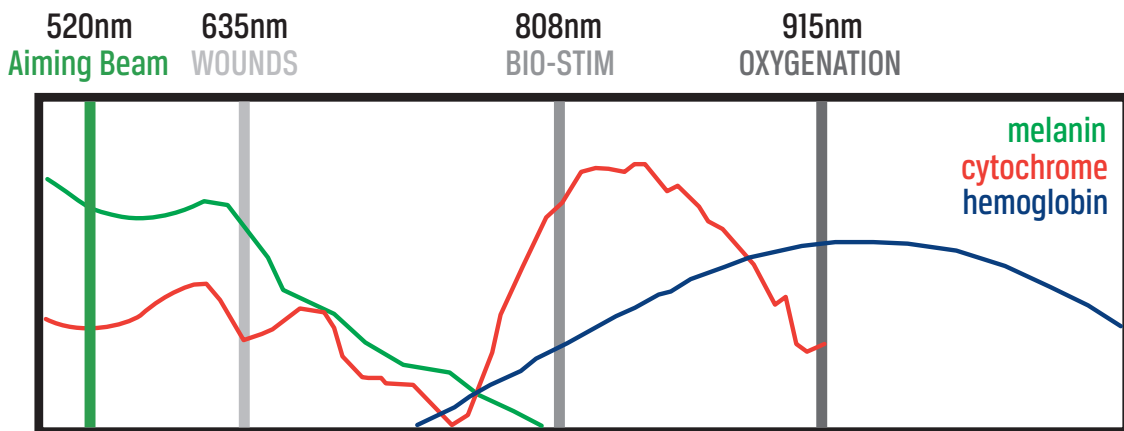

WOUNDS

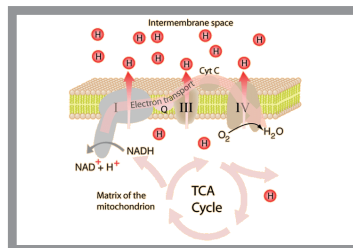

BIO-STIMULATION

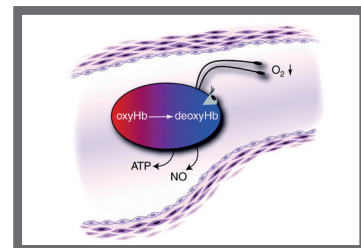

OXYGENATION

Infrared wavelengths are ideal for biostimulation. No study exists identifying superiority of one wavelength over another. By combining 635 nm, 808 nm and 915 nm ($\pm 10\%$), **doctorVet Plus** ensures optimal penetration and stimulation of targeted tissues without heating surface tissue. This allows **doctorVet Plus** to meet the broad treatment needs of veterinary practice.

MORE THAN LASER THERAPY.

- WiFi for Software Update, Client Communication, MWL and PACS Integration
- Patient Treatment Report
- doctorVet Metrics (DVM)
- Treatment Series Scheduling
- Treatment Reminders
- DICOM Medical Records
- PACS Integration
- Cloud-Based Data
- Modality Work List Integration
- PDT (Photodynamic Therapy)



Pain Management



Wound Healing



Biostimulation



Acupuncture



Oral Conditions



Dermatology

TOBI COMES WHEN YOU CALL.

Proper training is essential to effective clinical use. But it is difficult to prepare a nurse to treat every condition presented. **doctorVet Plus** teaches in real-time with step-by-step, on board instructions. Our team of experts, with a combined 20+ years of veterinary laser therapy experience, know how to get clinical results. Within every phase of treatment, specialized advice is displayed to guide the technician through the intricacies of treatment technique that provide maximum therapeutic efficacy.

TOBI (Treatment On Board Instructions) is a tutorial that guides you through each protocol phase to ensure success for all your doctorVet Plus treatments. The "i" instruction screen provides additional treatment advice to resulting in optimal treatment efficacy and training.



20 YEARS OF EXPERIENCE IN EVERY doctorVet Plus



80% of the patients that walk through your doors might benefit from doctorVet Plus treatments. doctorVet Plus provides a solution with the latest technology at the most affordable cost.

doctorVet Plus does this by incorporating features that have scientifically proven benefit, and omitting the high-cost features other companies try to convince you are needed. Focus on what the science says; Power, Pulsing, and Dosimetry are the only proven parameters.

Feature	doctorVet Plus	doctorVet V2
Peak Power	16 Watts	16 Watts
Wavelength Power	635 nm 500 mW 810 nm 8 W 915 nm 8 W	635 nm 500 mW 810 nm 8 W 915 nm 8 W
Color Touch Screen	●	●
Battery Power	●	●
5-Phase Treatment Protocols	Canine, Feline, Equine, Small Mammal, Avian, Exotic	Canine, Feline, Equine, Small Mammal, Avian, Exotic
TOBI	Multimedia	Multimedia
Massager TX Head	Included	Included
Spotter TX Head	Optional	Optional
Sweeper TX Head	Optional	Optional
Metrics Analysis	●	●
WiFi	●	●
Software Update	WiFi (automatic)	WiFi (automatic)
PDT (Photodynamic Therapy)	●	●
Patient Database	●	
Treatment Scheduling	●	
Client Communications	●	
PACS/MWL Connectivity	●	



Certified Quality System



UNI EN ISO 9001:2015
UNI EN ISO 13485:2012



LAMBDA SPA

Via dell'Impresa, 1
36040 BRENDOLA (VI) Italy
T +39 0444.349165
F +39 0444.349954
info@lambdaspa.com
lambdaspa.com